

nano-tRNAseq Kit

Our tRNA Sequencing Solution: the perfect blend of physical and digital methodologies with our nano-tRNAseq kit, a complete solution designed for streamlined tRNA multiplexed sequencing.

Product	Catalog no	Rxns.
nano-tRNAseq Kit with demultiplexing	NTRSQ-12	12

Shipping: Blue Ice and Dry ice

Storage Conditions: store components according to this manual

Shelf Life: 12 months

<u>Description</u>: nano-tRNAseq enables single molecule sequencing of native full-length tRNA using the Oxford Nanopore platform.

- (i) quantification of tRNA abundances,
- (ii) infer on tRNA modifications sites and
- (iii) detect modification "circuits".

Updated for RNA004 kit/RNA flow cell chemistry (2024).

Suitable for: Eukaryotic/prokaryotic cells and tissues with annotated genome.

What's Included:

- 1. **Library Prep Kit** Comprehensive tools for preparing your tRNA library ready for sequencing with Oxford Nanopore reagents.
- 2. **IT-Based Demultiplexing Tool** A powerful software tool that enables demultiplexing of your sequencing run, accessible for 40 days post-activation.

For Internal Research and Service Use Only. Not Intended for Diagnostic or Therapeutic Use.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK	Pag 1 of pag 21
---	-----------------

Table of Contents

4
5
7
.0

Kit storage information

	Quantity	Storage
-20°C components	1 box	-20°C
-80°C components	1 box	-80°C
USB pen drive with instructions on how to perform the demultiplexing step after sequencing.	1	RT

Additionally Required Material:

- RNA Clean & Concentrator™-5 (Zymo, cat. no. R1015 or R1016)
- Nuclease free water (NFW)
- Ethanol 95-100%
- RNase Inhibitor, Murine 40000 units/ml (New England Biolabs, cat. no. M0314L)
- o Agencourt AMPure XP for PCR Purification (Beckman Coulter, cat. no. A63881)
- QubitTM HS RNA Assay Kit (Thermo Fisher Scientific, Q10210)
- o QubitTM HS dsDNA Assay Kit (Thermo Fisher Scientific, Q32851/Q33230)
- Direct RNA Sequencing SQK-RNA004 (Oxford Nanopore)
- FLO-MIN004RA flow cell (Oxford Nanopore)
- o T4 DNA Ligase (New England Biolabs, cat. no. M0202M)
- NEBNext® Quick Ligation Reaction Buffer (New England Biolabs, cat. no. B6058S)

Optional Material:

- Agilent 2100 Bioanalyzer
- Agilent RNA 6000 Nano Kit (or equivalent) (Agilent Technologies, cat. no. 5067-1511)
- Agilent Small RNA Kit (Agilent Technologies, cat. no. 5067-1548)

INTRODUCTION

nano-tRNAseq is an Immagina's technology that marks a significant milestone in tRNA and translatomics research. It enables the comprehensive sequencing of tRNA molecules in their full-length, native state, allowing for the simultaneous assessment of tRNA abundances and modification status. The **nano-tRNAseq Kit** is the only product that allows for the sequencing of native full-length tRNA molecules with an easy, fast and robust workflow.

Studying tRNAs is particularly challenging due to their extensive post-transcriptional modifications and dynamic roles in translation. For these reasons, traditional sequencing methods such as sequencing by synthesis often fail to capture such complexities in tRNAs. Other NGS-based methods are complementary to nano-tRNAseq since those are cDNA and PCR based. Building upon the pioneering work of the Dr. Eva Novoa Lab at Centre for Genomic Regulation (CRG), Barcelona, our method bypasses the need for cDNA/PCR sequencing (https://doi.org/10.1038/s41587-023-01743-6)

The first step in nano-tRNAseq comprises small RNA enrichment from total RNA (1). After deacylation (2), the tRNAs are bound with adaptors and subsequently to barcodes (3), allowing for multiplexing of up to 6 samples. The native tRNA is then sequenced on the Oxford Nanopore Technologies (ONT) platform (4). As tRNA is transversing the pore, each base generates an electrical current, which is converted into the tRNA sequence. After data analysis (5), information on tRNA abundance, coverage and post-transcriptional chemical modifications are obtained (Figure 1).

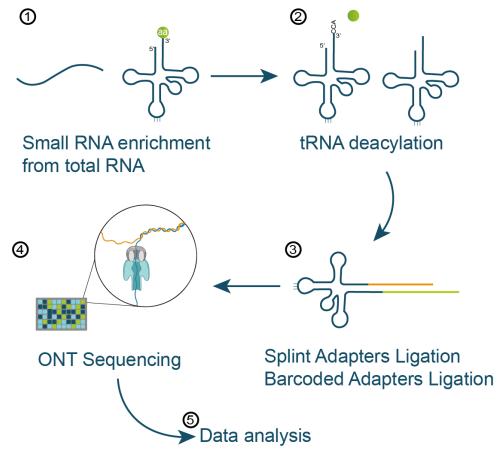


Figure 1. Simple overview of the nano-tRNAseq protocol.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK Pag	g 4 of pag 21
---	---------------

PROTOCOL WORKFLOW

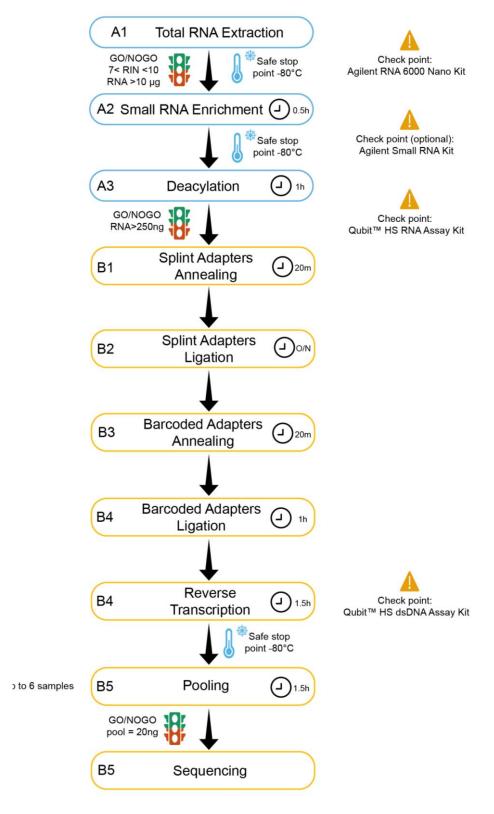


Figure 2. Overview of the nano-tRNAseq Kit workflow. In each box, steps are present in the left corner, while the time to complete the step is listed in the right corner (clock picture). Safe stopping points (thermometer pictures) and GO/NO-GO conditions (traffic lights) are placed right at the bottom of the corresponding steps. Optional and required checkpoints (caution signal) are also indicated on the right.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK	Pag 5 of pag 21
---	-----------------

Optimal Workflow Recommendations

- The nano-tRNAseq Kit has been optimized for the multiplexing of up to 6 samples per library.
- The multiplexing protocol allows you to process a total of 12 samples by generating pooled libraries composed of up to 6 samples each. If you wish, you can multiplex less samples per reaction. In case you want to multiplex less than 6 samples, we suggest you to use combinations of the barcoded adapters in the following order BC1 > BC2 > BC3 > BC4 > BC5 > BC6.
- Allocate at least 2 days for the completion of the entire workflow.
- Please make sure to purchase all the additionally required materials needed for the protocol before starting the experiment.
- This protocol is optimized for tRNA extraction starting from at least 10 μg of total RNA.
- At least 250 ng of small RNA are needed after the deacylation step to move forward with library preparation.
- Please be aware that the "Splinter Adaptor Ligation" reaction (step B2) can be performed starting from exactly 250 OR 500 ng of input.
- The final library pool must comprise exactly 20 ng of material, equally divided among the number of samples you decided to multiplex.
- The pooled tRNA library MUST be sequenced right away and cannot be stored for later processing.
- Before starting the sequencing, please select "Flow cell type": MIN004-RA, "Kit selection":
 Direct RNA Sequencing Kit, please deactivate "Basecalling" and be sure to select the saving
 of the .POD5 file formats.
- Note that at least 1.5M raw reads are needed to perform data analysis when multiplexing 6 samples.

A. tRNA EXTRACTION

Components and additional materials needed in this section:

Step N	Kit Component	Volume	Storage	Туре	Vial cap color
A3	Deacylation Buffer (DB)	1300 µL	-20°C	Vial	clear

Step N	Additional Material	Туре	
A1	Agilent RNA 6000 Nano Kit	Optional	
A2/A3	RNA Clean & Concentrator™-5	Required	
A1/2/3	Nuclease free water (NFW)	Required	
A1/2/3	Ethanol	Required	
A2	Agilent Small RNA Kit	Optional	
A3 Qubit™ HS RNA Assay Kit		Required	

Step A1. TOTAL RNA EXTRACTION

A1.1 Start extracting total RNA from your samples with any kit available at your facility. It is important to check for proper RNA integrity before proceeding with small RNA enrichment. The best approach is evaluating the extracted RNA by **Agilent 2100 Bioanalyzer** using the Agilent RNA 6000 Nano Kit (for non-eukaryotic species please use equivalent). An expected RNA Integrity Number (RIN) between 7 and 10 is needed for producing informative libraries.

NOTE: After total RNA extraction, it is important to enrich only the small RNA fraction to retrieve high quality tRNAs. Please proceed with **at least** 10 µg of total RNA to get enough material for library preparation.

Step A2. SMALL RNA (<200nt) ENRICHMENT

Small RNA Enrichment is achieved using the RNA Clean & Concentrator™-5 kit (Zymo catalog. no. R1015 or R1016). Perform all steps at room temperature and centrifugation at 12,000 g for 30 seconds, unless otherwise specified.

- A2.1 Add nuclease free water to the extracted RNA (at least 10 μg) from step A1.1 to a final total volume of 50 μL.
- □ A2.2 Prepare adjusted RNA Binding Buffer by mixing 50 μL of Zymo RNA Binding Buffer (ZBB) and 50 μL of ethanol (95-100%) for each sample you want to process.
- A2.3 Add 100 μL of adjusted RNA Binding Buffer (from step A2.2) to each sample and mix.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK Pag 7 of pag 2	Mod 080281	V1.2 D	Date Exp/rev 11/11/2024	Edited by MA IB MS EK		Pag 7 of pag 21
--	------------	--------	-------------------------	-----------------------	--	-----------------

- □ A2.4 Transfer the mixture to the Zymo-Spin™ Column and centrifuge. Save the flow-through: Small RNAs (17-200 nt) are in the flow-through!
- A2.5 Add 150 μL of ethanol and mix. Transfer the mixture to a new column and centrifuge. Discard the flow-through.
- A2.6 Add 400 μL Zymo RNA Prep Buffer to the column and centrifuge. Discard the flow-through.
- A2.7 Add 700 μL Zymo RNA Wash Buffer to the column and centrifuge. Discard the flow-through.
- A2.8 Add 400 μL Zymo RNA Wash Buffer to the column and centrifuge for 1 minute to ensure complete removal of the wash buffer. Carefully, transfer the column into a new RNase-free tube.
- A2.9 Add 11 μL of nuclease-free water directly to the column matrix, wait for 1 minute at RT and centrifuge.
- □ A2.10 The small RNAs are present in the flow-through. Keep the reaction tube containing the flow through.



SAFE STOPPING POINT (store at -80°C)

Note that, while the use of the RNA Clean & Concentrator[™]-5 kit (Zymo catalog. no. R1015 or R1016) is suggested, it is anyway possible to extract small RNA with any commercially available kit at your facility.

It is suggested to check the small RNA quality after extraction using a Bioanalyzer for small RNA (Agilent Small RNA kit).

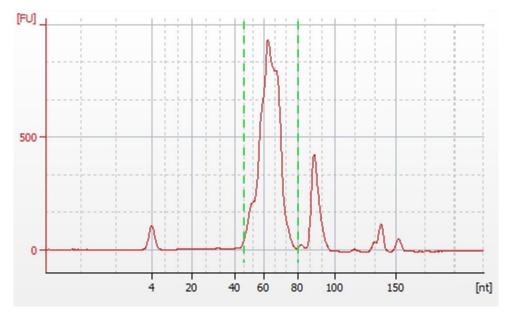


Figure 3. example of small RNA Bioanalyzer trace after RNA extraction and small RNA enrichment. A major peak between 50-80 nt is expected.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK Pag 8 of pag 21	
---	--

Step A3. DEACYLATION

	A3.1 Add 90 μL of deacylation buffer (DB) to the small RNA fraction from step A2.10 so that the total volume is 100 μL .
	A3.2 Incubate the reaction 30 minutes at 37°C.
R1	rify deacylated small RNA using the RNA Clean & Concentrator™-5 kit (Zymo catalog. no. 015 or R1016). Perform all steps at room temperature and centrifugation at 12,000 g for 30 conds, unless otherwise specified.
	A3.3 Add 200 µL of Zymo RNA Binding Buffer (ZBB) to each sample and mix.
	A3.4 Add 400 µL of ethanol 100% to each sample and mix well.
	A3.5 Transfer the mixture to the Zymo-Spin™ Column and centrifuge.
	Discard the flow-through.
	A3.6 Add 400 µL Zymo RNA Prep Buffer to the column and centrifuge.
	Discard the flow-through.
	A3.7 Add 700 µL Zymo RNA Wash Buffer to the column and centrifuge.
	Discard the flow-through.
	A3.8 Add 400 µL Zymo RNA Wash Buffer to the column and centrifuge for 1 minute to ensure complete removal of the wash buffer. Carefully, transfer the column into a new RNase-free tube.
	A3.9 Add 8 μL of nuclease-free water directly to the column matrix, wait for 1 minute and centrifuge.
	A3.10 The small RNAs are present in the flow-through.
	Keep the reaction tube containing the flow through.
	A3.11 Quantify 1 µL from A3.10 using a Qubit™ HS RNA Assay Kit. The recovered material should be at least 250 ng of deacylated small RNA.
	Please use deacylated small RNA for library preparation (step C) within the same day.
	Fricase use deacyrated sinan Kiva for horary preparation (step c) within the same day.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK Pag 9 of pag 21

B. tRNA LIBRARY PREPARATION

nano-tRNAseq Kit components needed in this section:

Step N	Kit component	Volume	Storage	Туре	Vial cap color
B1/B3	Annealing buffer 1 (AB1)	50 μL	-20°C	vial	Yellow
B1/B3	Annealing Buffer 2 (AB2)	50 μL	-20°C	vial	Yellow
B1	Splint adapter A (SA)	50 μL	-80°C	strip	
B1	Splint adapter B (SB)	50 μL	-80°C	strip	
B2	PEG 8000 (PEG)	300 µL	-20°C	vial	Yellow
B2	Buffer T1 (BT1)	50 μL	-20°C	vial	Yellow
B2	T1 enzyme (T1)	27 µL	-20°C	vial	Yellow
B3	Barcoded Adapter BC1 FWD (1 F)	25 µL	-80°C	strip	
B3	Barcoded Adapter BC1 REV (1 R)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC2 FWD (2 F)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC2 REV (2 R)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC3 FWD (3 F)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC3 REV (3 R)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC4 FWD (4 F)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC4 REV (4 R)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC5 FWD (5 F)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC5 REV (5 R)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC6 FWD (6 F)	25 μL	-80°C	strip	
B3	Barcoded Adapter BC6 REV (6 R)	25 μL	-80°C	strip	
B4	Buffer T2 (BT2)	65 µL	-20°C	vial	Red
B4	T2 enzyme (T2)	20 μL	-20°C	vial	Red
B4	dNTPs	30 µL	-20°C	vial	Green
B4	Buffer T3 (BT3)	110 µL	-20°C	vial	Green
B4	T3 enzyme (T3)	10 μL	-20°C	vial	Green

INIOQ USUZSI	Mod 080281	V1.2	Date Exp/rev 11/11/2024	Edited by MA IB MS EK		Pag 10 of pag 21
--------------	------------	------	-------------------------	-----------------------	--	------------------

Additional materials needed in this section:

Step N	Additional Material	Туре
B1/B2/B3/B4	RNase Inhibitor, Murine	Required
B2/B4	Agencourt AMPure XP	Required
B2/B4	Ethanol	Required
B2/B3/B4/B5	Nuclease Free Water (NFW)	Required
B4	Qubit [™] HS dsDNA Assay Kit	Required
B5	Direct RNA sequencing SQK-RNA004 (Oxford Nanopore)	Required
B5	FLO-MIN004RA flow cell (Oxford Nanopore)	Required

Step B1. SPLINT ADAPTERS ANNEALING

Please note that a single 10 μ L reaction for the annealing of splint adapters (B1.1) will be sufficient to handle 5 or 2 samples according to the input RNA that you wish to use in reaction B2.1 (250 or 500 ng, respectively), so you may need to perform more than one reaction from section B1.1 depending on the number of samples you decide to process in a single experiment.

□ B1.1 Mix the following reagents in a 1.5 mL reaction tube. Please note that the volumes in Table 1 are intended for one single reaction of splint adapters annealing. Plan the number of reactions according to the number of samples that you wish to process, and the RNA input you plan to use in reaction B2.1.

Reagent	Volume (µL)
Annealing Buffer 1 (AB1)	1
Annealing Buffer 2 (AB2)	1
Splint Adapter A (SA)	3.75
Splint Adapter B (SB)	3.75
RNase Inhibitor, Murine	0.5
Total volume	10

Table 1. Components' volumes to use in step B1 in one single reaction of splint adapters annealing.

□ B1.2 Mix the reactions well by pipetting, then heat to 75°C for 15 sec, and ramp down to 25°C at 0.1°C/s. Store on ice until further use. Once used, please toss the leftovers.

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK Pag 12 of pag 21
--

Step B2. SPLINT ADAPTERS LIGATION

NOTE: Start the library preparation with 250 or 500 ng of deacylated small RNA. Adjust the amount of annealed Splint Adapter to use according to the input material (see Table 2 below).

B2.1 Mix the following reagents in a 1.5 mL reaction tube. For clarity, volumes indicated in Table 2 are to be considered for one reaction only and must be repeated for each sample.

	250 ng of RNA	500 ng of RNA
Reagent	Amount (µL)	Amount (µL)
Deacylated RNA from step A3.10	<mark>250 ng</mark> (Χ μL)	<mark>500 ng</mark> (Χ μL)
Annealed Splint Adapter from step B1.2	1.9	3.8
PEG 8000	10	10
Buffer T1 (BT)	2.5	2.5
T1 Enzyme (T1)	2	2
RNase Inhibitor, Murine	0.5	0.5
H ₂ O	8.1 - X	6.2 - X
Total volume	25	25

Table 2. Components' volumes to use in step B2 for one reaction. The two columns show the volumes of reagents according to the input of small RNA deacylated chosen (250 ng or 500 ng).

addording to the input of small NW deady lated chosen (200 mg of 600 mg).
B2.2 Incubate the reaction overnight at 4°C.
B2.3 Purify the reaction by adding 50 μ L of vortexed Agencourt AMPure XP beads (warm the beads at RT for 30 minutes before use) to each sample and mix well by pipetting the entire volume up and down at least 10 times.
B2.4 Incubate at RT for 10 minutes.
B2.5 Place the tubes on a magnetic rack and discard the supernatant when clear and colorless
tRNAs are now attached to the beads!
B2.6 Keep the tubes on the magnetic rack. Add 200 μL of EtOH 70% freshly prepared to the beads. Incubate for 30 seconds and then remove the supernatant.
beads. Incubate for 30 seconds and then remove the supernatant.

Avoid over-drying the beads (pellet cracked) as this will significantly decrease elution efficiency.

- B2.9 Remove the tubes from the magnetic rack and resuspend the beads in 9 μL of nuclease-free water. Mix thoroughly by pipetting up and down to ensure complete bead dispersion. Incubate at RT for 10 minutes.
- □ B2.10 Place the tubes on a magnetic rack until the solution is completely clear.

tRNAs are now in the supernatant!

B2.11 Remove the 9 μL of supernatant and place into a clean 0.2 mL nuclease-free tube.

Mod 080281	V1.2	Date Exp/rev 11/11/2024	Edited by MA IB MS EK	Pag 13 of pag 21
11100 000201	* +	Date Exp, (et 11, 11, 202)	Laitea by itin tib itis Lit	1 ag 13 01 pag 21

Step B3. BARCODED ADAPTERS ANNEALING

NOTE: Each pair of barcoded adapters (FWD and REV) needs to be annealed following the passages below. For clarity, Table 3 shows an example for the barcode pair **1** F and **1** R. Please, repeat the reaction for every pair of barcodes needed (up to 6 per library preparation).

□ B3.1 Mix the following reagents in a 1.5 mL reaction tube:

Reagent	Volume (µL)
Annealing Buffer 1 (AB1)	1
Annealing Buffer 2 (AB2)	1
Barcoded Adapter BC1 FWD (1 F)	1.5
Barcoded Adapter BC1 REV (1 R)	1.5
H ₂ O	4.5
RNase Inhibitor, Murine	0.5
Total volume	10

Table 3. Components' volumes to use in step B3.

□ B3.2 Mix the reactions well by pipetting, then heat to 75°C for 15 sec, and ramp down to 25°C at 0.1°C/s. Store on ice until further use. Once used, please toss the leftovers.

Step B4. BARCODED ADAPTERS LIGATION AND REVERSE TRANSCRIPTION

□ B4.1 Mix the following reagents in a 0.2 mL nuclease-free reaction tube. Please perform separate reactions for each sample/barcoded adapter you are processing.

Reagent	Volume (µL)
Buffer T2 (BT2)	3
tRNA from step B2.11	8.5
Annealed Barcoded Adapter* from step B3.2	1.5
RNase Inhibitor, Murine	0.5
T2 Enzyme (T2)	1.5
Total volume	15

^{*}Please use different Barcoded Adapters for different samples

Table 4. Components' volumes to use in step B4.1. Please perform separate reactions for each sample/barcoded adapters you are processing.

- □ B4.2 Incubate the reaction for 10 minutes at RT.
- B4.3 Meanwhile, prepare the reverse transcription master mix as follows. Please consider that the volumes in Table 5 are to be considered for one single sample and corresponding barcoded adapters.

Reagent	Volume (μL)
H ₂ O	14.5
Buffer T3 (BT3)	8
dNTPs	2
Total volume	24.5

Table 5. Components' volumes to use in step B4.3. Volumes in this table are to be considered for one single sample and corresponding barcoded adapters.

- □ B4.4 Add the master mix to the reaction tube containing the barcoded adapters-ligated tRNA from step B4.2. Mix by pipetting.
- B4.5 Add 0.5 μL of T3 enzyme (T3) to the reaction and mix by pipetting.
- □ B4.6 Incubate at 60°C for 30 minutes, then at 85°C for 5 minutes, and bring to 4°C before proceeding with the next step.
- B4.7 Purify the reaction by adding 80 μL of vortexed Agencourt AMPure XP beads (warm the beads at RT for 30 min before use) to each sample and mix well by pipetting the entire volume up and down at least 10 times.
- □ B4.8 Incubate at RT for 10 minutes.
- B4.9 Place the tubes on a magnetic rack and discard the supernatant when clear and colorless.

NA 1.000204		D : 5 / 44/44/2024	5 I'' II AAA IB AAC 514	D 45 (24
Mod 080281	V1.2	Date Exp/rev 11/11/2024	Edited by MA IB MS EK	Pag 15 of pag 21

tRNAs are now attached to the beads!

B4.10 Keep the tubes on the magnetic rack. Add 200 μ L of EtOH 70% freshly prepared to the beads. Incubate for 30 seconds and then remove the supernatant.
B4.11 Repeat the washing step.
B4.12 Remove supernatant and let the beads pellet dry on the magnetic rack at room temperature for ~2-4 minutes. Avoid over-drying the beads (pellet cracked) as this will significantly decrease elution efficiency.
B4.13 Remove the tubes from the magnetic rack and resuspend the beads in 8.5 μ L of nuclease-free water. Mix thoroughly by pipetting up and down to ensure complete bead dispersion. Incubate at room temperature for 10 minutes.
B4.14 Place the tubes on a magnetic rack until the solution is completely clear.
tRNAs are now in the supernatant!
B4.15 Remove the 8.5 μL of supernatant and place into a clean 1.5 mL reaction tube.
B4.16 Quantify 2 μL of the eluate using a Qubit TM HS dsDNA Assay Kit.

SAFE STOPPING POINT (store at -80°C for up to one week)

Step B5. RNA POOLING, RMX LIGATION AND SEQUENCING

NOTE: At this step it is possible to decide how many samples to pool together as long as 20 ng of total material is reached.

B5.1 Pool in a 1.5 mL tube the barcoded samples from step B4.16 so that the total RNA amount is 20 ng in 23 μL (if necessary, use nuclease-free water to reach the requested volume). Use the following table as a guide:

Samples barcoded with	Quantity (ng)	Volume (μL)
BC1	3.3	
BC2	3.3	
BC3	3.3	
BC4	3.3	
BC5	3.3	
BC6	3.3	
NFW	if needed	
Total	20	23

□ C5.2 Proceed with the <u>Library preparation from Oxford Nanopore SQK-RNA004</u> starting from step 24 of section 3 ("Library Preparation").


⚠ IMPORTANT: The tRNA library must be sequenced immediately and cannot be stored for later use.

MIMPORTANT: BEFORE STARTING THE SEQUENCING

⚠ Select "Flow cell type": FLO-MIN004-RA



Select "Kit selection": Direct RNA Sequencing Kit (SQK-RNA004).

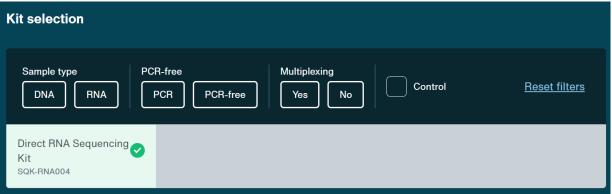

Deactivate "Basecalling".



A Select the saving of the .POD5 file format in the software before starting the sequencing.



⚠ All the information related to the demultiplexing step are present in the USB pen drive that has been provided with the kit. Please, make sure to read all the instructions before loading the library on the flowcell.

Contacts



General information

info@immaginabiotech.com

Sales support (quoting, ordering, and order status update)

orders@immaginabiotech.com

Technical service (technical inquiries and quality complaints)

techsupport@immaginabiotech.com



Viale dell'industria, 47, 38057, Pergine Valsugana (TN), ITALY



www.immaginabiotech.com



+39 04611787270

Immagina's mission is to develop unique and smart enabling technologies to break down the walls in translational studies. Please visit our website https://immaginabiotech.com/ for a complete overview of our products & services and our proprietary technologies.

Notes:			

Mod 080281 V1.2 Date Exp/rev 11/11/2024 Edited by MA IB MS EK	Pag 20 of pag 21
---	------------------

nano-tRNAseq Kit	

Mod 080281	V1.2	Date Exp/rev 11/11/2024	Edited by MA IB MS EK	Pag 21 of pag 21
10100 000201	V 1.2	Date Exp/164 11/11/2024	Luited by WIA ID WIS ER	1 46 21 01 946 21